

Eukaryotic Cluster Of Differentiation 55 (CD55)

Catalog No.: TD-TP04848 50µg

Sequence Information

Species: Human

Gene ID:1604

Swiss Prot:P08174

Synonyms:DAF; CR; TC; Complement

decay-accelerating factor; Decay

Accelerating Factor

Residues:Asp35-Ser353

DCGLPPDVPNAQPALEGRSTFPEDTVITYKCEESFVKIPGEKDSVICKGSQWS

DIEEFCNRSCEVPTRLNSASLKQPYITQNYFPVGTVEYECRPGYRREPSLSPK

LTCLQNLKWSTAVEFCKKKSCPNPGEIRNGQIDVPGGILFGATISFSCNTGYKL

FGSTSSFCLISGSSVQWSDPLPECREIYCPAPPQIDNGIIQGERDHYGYRQSVT

YACNKGFTMIGEHSIYCTVNNDEGEWSGPPPECRGSLTSKVPPTVQKPTTVNV

PTTEVSPTSQKTTTKTTTPNAQATRSTPVSRTTKHFHETTPNKGSGTTS

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: N-terminal His-Tag

Subcellular Location: Membrane, Secreted.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 150µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.2

Predicted Molecular Mass: 36.6kDa

Accurate Molecular Mass: 58kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

电话: 027-87561633 手机: 189-8623-2100

QQ客服: 121545059 网址: <http://www.itiande-neuron.com>

公司地址: 武汉东湖新技术开发区高新大道666号光谷生物城生物创新园B1栋b274-9室

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

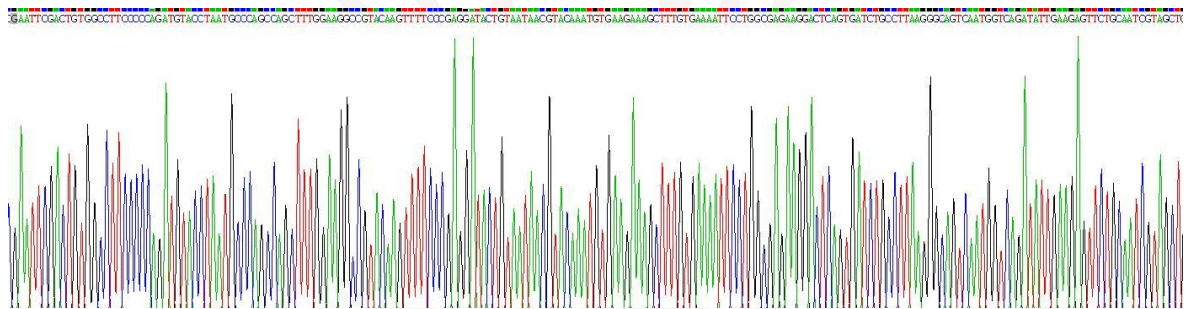


Figure 1. Gene Sequencing (Extract)

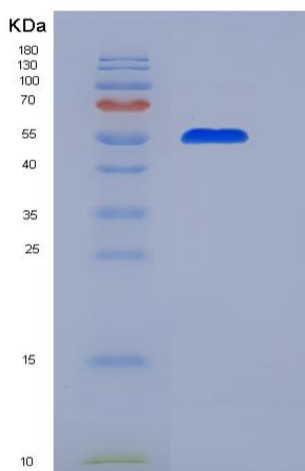


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.